

Advance Plus

Touch Free Intermittent Catheter System


Confidence. Pure and Simple.



Tear here to empty

Hollister Advance Plus Intermittent Catheter System features *Touch Free* technology enabling users to catheterize with confidence — anytime, anywhere — supporting user mobility and independence. The unique *Touch Free* protective introducer tip helps reduce the risk of carrying bacteria up into the urinary system, a cause of urinary tract infections. An integrated bag shields the sterile catheter from environmental contamination before and during insertion.

 Hollister

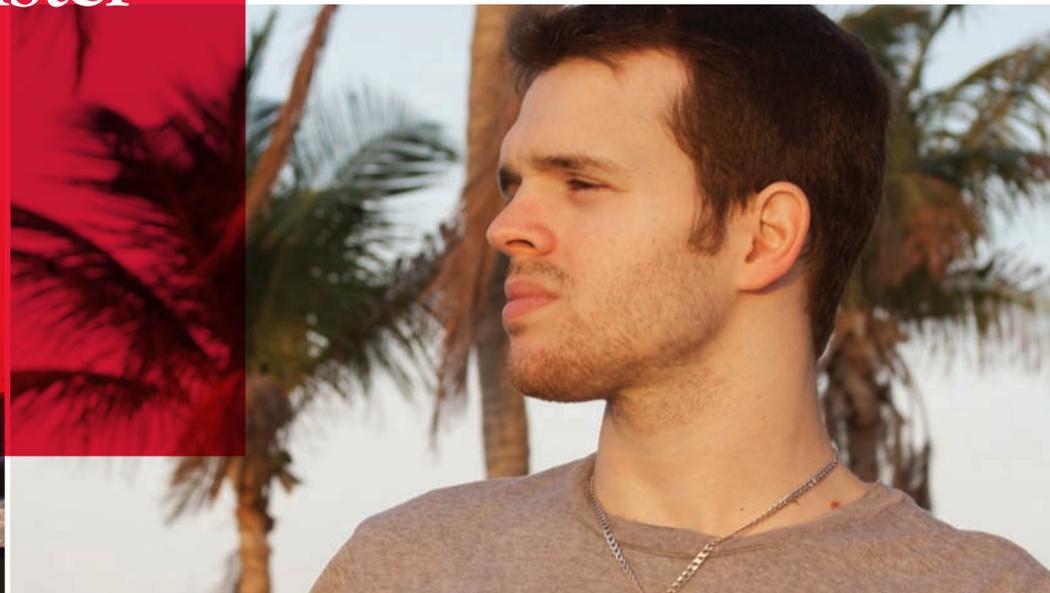

Advance Plus

Touch Free Intermittent Catheter System

1 Touch Free protective introducer tip

- protects catheter tip from exposure to initial 15 mm of the urethral opening
- helps reduce the risk of carrying bacteria up into the urinary system
- serves as a guide for proper and easy insertion

2 Touch Free 1500 mL collection bag

- shields catheter from possible contamination before and during insertion
- convenient bag handle makes for easy handling, even for users with limited manual dexterity
- ensures secure urine collection
- supports user independence—catheterize anytime, anywhere

3 Ultra-smooth catheter eyelets

- helps to ensure gentle insertion and withdrawal
- enhances user comfort

4 Patented user-regulated gel reservoir

- no additional lubrication required
- enhances user convenience and comfort
- stabilizes catheter for easy insertion and withdrawal

5 Red catheter ring cap


- helps keep protective introducer tip clean and protected after package is open
- finger hole makes it easy to remove and replace, even with limited dexterity
- red oval shape doesn't roll if dropped

Coudé and Straight insertion tips

- accommodates user preference and need

Latex- and phthalates-free PVC material

- safe for users with latex-related allergies
- environmentally responsible

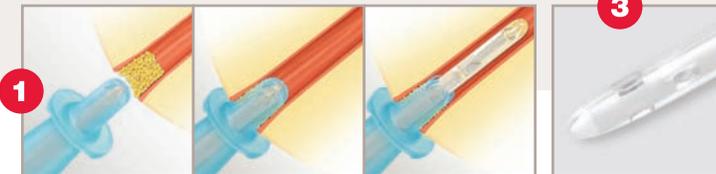

A Hollister Technology

Advance Plus Touch Free Intermittent Catheter System

Stock No	Tip	Size (16 in, 40cm)	Box Qty
94064	Straight	06Fr	100
94084	Straight	08Fr	100
94104	Straight	10Fr	100
94124	Straight	12Fr	100
94144	Straight	14Fr	100
94164	Straight	16Fr	100
94184	Straight	18Fr	100
95124	Coudé	12Fr	100
95144	Coudé	14Fr	100
95164	Coudé	16Fr	100

Advance Plus Kit Touch Free Intermittent Catheter

Stock No	Tip	Size (16 in, 40cm)	Box Qty
96064	Straight	06Fr	100
96084	Straight	08Fr	100
96104	Straight	10Fr	100
96124	Straight	12Fr	100
96144	Straight	14Fr	100
96164	Straight	16Fr	100
96184	Straight	18Fr	100
97124	Coudé	12Fr	100
97144	Coudé	14Fr	100
97164	Coudé	16Fr	100

Kit includes: Advance Plus Intermittent Catheter, gloves, underpad drape, and Benzalkonium swabsticks



Warning: To help reduce the potential risk of infection and/or other complications, do not re-use. Dispose of appropriately after procedure. If discomfort or any sign of trauma occurs, discontinue use immediately and consult your healthcare professional.

Caution: Federal (USA) law restricts this device to sale by or on the order of a licensed healthcare professional. Prior to use of this device, be sure to read (i) the complete information on how to use this device including Warnings, Cautions, and Instructions for Use, and (ii) all other package inserts and labels supplied with the product and accessories.

Hollister and logo, and Advance Plus are trademarks of Hollister Incorporated.
©2010 Hollister Incorporated. 911158-710

Please consult a medical professional before using this product if any of the following conditions are present: Severed urethra, unexplained urethral bleeding, pronounced stricture, false passage, urethritis – inflammation of the urethra, prostatitis – inflammation of the prostate gland, epididymitis – inflammation of the epididymis (testicle tube). Self-catheterization should only be carried out after medical advice and only in accordance with the instructions provided. Always follow the care plan and advice given by your healthcare professional. For urethral intermittent self-catheterization (ISC), it is typical to catheterize at least 4 times a day at intervals of 6 to 8 hours. If you are unsure about your catheterization, please contact your regular healthcare professional.

Hollister Incorporated
2000 Hollister Drive
Libertyville, Illinois 60048
1.888.740.8999

Distributed in Canada by
Hollister Limited
95 Mary Street
Aurora, Ontario L4G 1G3
1.800.263.7400

www.hollister.com